Shower benches

All our shower benches are ideal for usage in the shower, bathtub, bathroom or other humid environments. The specially designed rubbers at the end of each leg will create a safe stability for the bench and avoid the sliding of the bench.

All models have leg adjustments. You can easily adjust the length of the legs by pressing the snap buttons and sliding the tube in- or outwards. Make sure the snap button securely locks into the hole. Ensure that all four legs are of the same height before using your shower bench.

Ask assistance if you are uncertain on the usage of your product.

Users weight: 150 kg

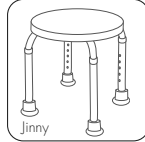
The shower benches can be cleaned with normal household cleaning products. Avoid using aggressive cleaning products which can cause discoloration.

Dimensions Jill

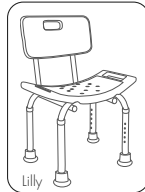
width: 51 cm
length: 42 cm
seat height: 35 - 53 cm
seat width: 51 cm
leg adjustment: 15 cm
weight: 2 kg

**Dimensions Jinny**

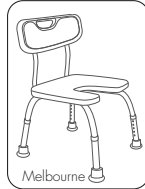
width: 35 cm
length: 35 cm
seat height: 33 - 51 cm
seat width: 32 cm
leg adjustment: 15 cm
weight: 1,5 kg

**Dimensions Lilly**

width: 51 cm
length: 47 cm
height: 77 - 87 cm
height backrest: 38 cm
seat width: 51 cm
seat height: 36 - 46 cm
leg adjustment: 15 cm
weight: 3 kg

**Dimensions Melbourne**

width: 48 cm
length: 56 cm
height: 72 - 82 cm
height backrest: 37 cm
seat width: 41 cm
seat height: 34 - 44 cm
leg adjustment: 15 cm
weight: 2,8 kg

**Dimensions Sidney**

width: 50 cm
length: 52 cm
seat height: 34 - 44 cm
seat width: 41 cm
leg adjustment: 15 cm
weight: 2 kg

Subject to technical alterations. © 2016

VERMEIREN

Opening up new horizons in mobility

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